


NOVA

BEST LENS FOR OFFICE WORK

OFFICE

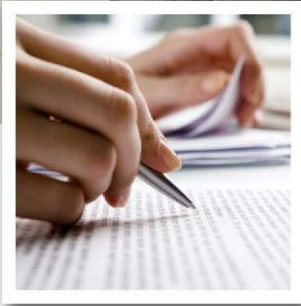

The best lens for office work

NOVA OPTICAL LABS
TRULY ALBERTAN INDEPENDENT

CALGARY 1-403-272-2007

RED DEER 1-403-346-0999

LETHBRIDGE 1-403-329-0041

WWW.NOVAOPTICALLABS.CA

NOVA Office freeform

When Your Glasses Are Not Enough

Lenses with only a reading area provide clear vision of objects located close up, while progressive lenses only provide a narrow intermediate vision.


Office Reader is an ideal solution for near and intermediate vision, solving most limitations of other glasses. Ideal for presbyopes with special needs for intermediate and near vision such as office workers, cooks, artists and many other occupations.



The occupational lens exclusive for both intermediate and near vision

NOVA OFFICE design overview


Nova Office is the best occupational lens. It offers wide near and intermediate visual areas combined with a minimum lateral astigmatism. This perfect union reduces the adaptation period to almost immediate; providing the wearer with great comfort and allowing a natural posture when working on computers. Several design choices mean that you can prescribe the right option for your patient no matter what the wear's prescription or visual needs.

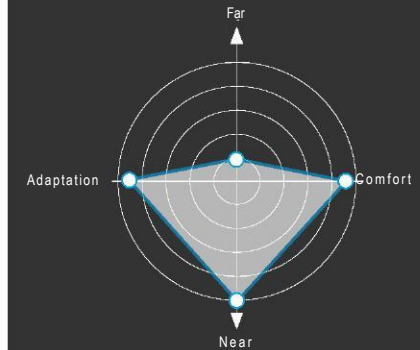

Select the option that fits with your needs

Nova Office offers several clear vision depths that provide the patients with a visual solution perfectly adapted to their individual needs. Available in three different

options; the user will experience a new visual sensation due to the wide intermediate and near vision combined with exceptional comfort.

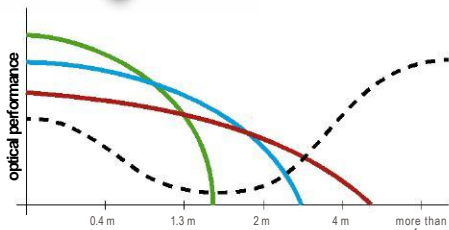

Relationship between vision performance and distance to the object

■ Reader 1.3 m	Up to 1.3 meters (4 ft) of clear vision
■ Reader 2 m	Up to 2 meters (6.5 ft) of clear vision
■ Reader 4 m	Up to 4 meters (13 ft) of clear vision
■ Progressive	variable performance



When to use my Nova Office lenses?

Nova Office Lens should be used when having special Requirement for near and intermediate vision.

For normal distance activities such as driving, sports, and everyday wear other design choices are preferred.



Enhanced Near



Enhanced For Computers



Not Suitable For Driving

NOVA OFFICE TECHNICAL DATA

Type	Occupational lens		
Performance	Offers unparalleled performance compared with other lenses intermediate and near visual fields.		
Vision Depth	Adaptable		
3 Available Options	Reader 1.3 m – clear vision from 33 cm (12 in) up to 1.3m (4 FT) Reader 2 m – clear vision from 33 cm (12 in) up to 2 m (6.5FT) Reader 4 m – clear vision from 33 cm (12 in) up to 4 m (13 ft)		
Minimum Fitting Height	18 mm	(14MM CORRIDOR)	2MM DROP
Note	not recommended for use while driving.		

- AVAILABLE IN;
- 1.50-CLEAR
 - 1.53- TRIVEX CLEAR
 - 1.58 - QUATREX CLEAR
 - 1.60 – CLEAR
 - 1.67- CLEAR
 - 1.70- CLEAR
 - 1.74- CLEAR

